

SATB2 (DNA-Binding Protein SATB2) (EP281)

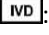
Rabbit anti-human SATB2 (DNA-Binding Protein SATB2) Monoclonal Antibody (Clone EP281)

REFERENCES AND PRESENTATIONS¹

- **ready-to-use (manual or LabVision AutoStainer)**
MAD-000747QD-3
MAD-000747QD-7
MAD-000747QD-12
- **Ready-to-use (MD-Stainer)²**
MAD-000747QD-3/V
MAD-000747QD/V
- **concentrated**
MAD-000747Q - 1:50 recommended dilution

COMPOSITION

Anti-human SATB2 rabbit monoclonal antibody purified from serum and prepared in 10mM PBS, pH 7.4, with 0.2% BSA and 0.09% sodium azide

INTENDED USE : Immunohistochemistry (IHC) on paraffin embedded tissues. Not tested on frozen tissues or Western-Blotting

CLONE: EP281³

Ig ISOTYPE: Rabbit IgG

IMMUNOGEN: A synthetic peptide corresponding to residues of human SATB2 protein.

SPECIES REACTIVITY: In vitro diagnostics in humans. Not tested in other species

DESCRIPTION AND APPLICATIONS:

DNA-binding protein SATB2, also known as Special AT-rich sequence-binding protein 2, is a nuclear matrix-associated transcription factor. SATB2 acts as a docking site for chromatin remodelling enzymes and recruits co-activators and co-repressors to control nuclear gene expression. SATB2 also regulates skeletal development, osteoblast differentiation, and modulates immunoglobulin expression.

In normal tissues, strong nuclear SATB2 expression is observed in essentially all glandular cells lining the lower gastrointestinal tract, including the appendix,



colon, and rectum. SATB2 is also expressed in a subset of neuronal cells from the cerebral cortex and hippocampus. In tumor tissues, SATB2 is detected in cancer cells of colorectal origin and may function as a clinically useful diagnostic marker for colorectal cancer (CRC). In a multi-cohort study with 1882 primary and metastatic CRCs, SATB2 shows high sensitivity (85%) for CRC, and further enhanced to 93% when stained in conjunction with Cytokeratin 20. A recent study showed SATB2 expression in 89% of medullary carcinomas of the large intestine. SATB2 has been suggested as a valuable prognostic marker: high SATB2 expression was determined as an independent marker of good prognosis and sensitivity to chemotherapy and radiation in CRC while loss of SATB2 expression was correlated with poor prognosis in laryngeal carcinoma patients.

IHC POSITIVE CONTROL: Colon for normal tissue and colorectal cancer for abnormal tissue.

VISUALIZATION: Cell nuclei

IHC RECOMMENDED PROCEDURE:

- 4µm thick section should be taken on charged slides; dry overnight at 60°C
- Deparaffinise, rehydrate and HIER (heat induced epitope retrieval) – boil tissue in the Pt Module using Vitro S.A EDTA buffer pH8⁴ for 20 min at 95°C. Upon completion rinse with 3-5 changes of distilled or deionised water followed by cooling at RT for 20 min
- Endogenous peroxidase block - Blocking for 10 minutes at room temperature using peroxidase solution (ref. MAD-021540Q-125)
- Primary antibody: incubate for 10 minutes [The antibody dilution (when concentrated) and protocol may vary depending on the specimen preparation and specific application. Optimal conditions should be determined by the individual laboratory]
- For detection use Master Polymer Plus Detection System (HRP) (DAB included; ref. MAD-000237QK)
- Counterstaining with haematoxylin and final mounting of the slide

STORAGE AND STABILITY:  Stored at 2-8°C. Do not freeze.  Once the packaging has been opened it can be stored until the expiration date of the reagent indicated on the label. If the reagent has been stored

¹ These references are for presentation in vials of Low Density Polyethylene (LDPE) dropper. In case the products are used in automated stainers, a special reference is assigned as follows:

- / L: Cylindrical screw-cap vials (QD-3 / L, QD-7 / L, QD-12 / L).

- / N: Polygonal screw-cap vials (QD-3 / N, QD-7 / N, QD-12 / N).

For different presentations (references / volumes) please contact the supplier.

² For Technical specifications for MD-Stainer, please contact your distributor.

³ SATB2 clone EP281 is manufactured using Epiteomics's RabMab® technology under U.S. Patent Nos. 5,675,063 and 7,402,409

⁴ Ref: MAD-004072R/D



under other conditions to those indicated in this document, the user must first check its correct performance taking into account the product warranty is no longer valid.

WARNINGS AND PRECAUTIONS:

1. Avoid contact of reagents with eyes and mucous membranes. If reagents come into contact with sensitive areas, wash with copious amounts of water.
2. This product is harmful if swallowed.
3. Consult local or state authorities with regard to recommended method of disposal.
4. Avoid microbial contamination of reagents.

SAFETY RECOMMENDATIONS

This product is intended for laboratory professional use only. The product is NOT intended to be used as a drug or for domestic purposes. The current version of the Safety Data Sheet for this product can be downloaded by searching the reference number at www.vitro.bio or can be requested at regulatory@vitro.bio.

BIBLIOGRAPHY

1: Berg KB, Schaeffer DF. SATB2 as an Immunohistochemical Marker for Colorectal Adenocarcinoma: A Concise Review of Benefits and Pitfalls. Arch Pathol Lab Med. 2017 Oct;141(10):1428-1433.

2: Yang C, Sun L, Zhang L, Zhou L, Zhao M, Peng Y, Niu D, Li Z, Huang X, Kang Q, Jia L, Lai J, Cao D. Diagnostic Utility of SATB2 in Metastatic Krukenberg Tumors of the Ovary: An Immunohistochemical Study of 70 Cases With Comparison to CDX2, CK7, CK20, Chromogranin, and Synaptophysin. Am J Surg Pathol. 2017 Sep 13. doi: 10.1097/PAS.0000000000000951. [Epub ahead of print]

3: Giannico GA, Gown AM, Epstein JI, Revetta F, Bishop JA. Role of SATB2 in distinguishing the site of origin in glandular lesions of the bladder/urinary tract. Hum Pathol. 2017 Sep;67:152-159

4: Matsushima J, Yazawa T, Suzuki M, Takahashi Y, Ota S, Nakajima T, Yoshino I, Yokose T, Inoue T, Kawahara K, Nakatani Y. Clinicopathological, immunohistochemical, and mutational analyses of pulmonary enteric adenocarcinoma: usefulness of SATB2 and β -catenin immunostaining for differentiation from metastatic colorectal carcinoma. Hum Pathol. 2017 Jun;64:179-185.

5: Li Z, Roth R, Rock JB, Lehman A, Marsh WL, Suarez A, Frankel WL. Dual Immunostain With SATB2 and CK20 Differentiates Appendiceal Mucinous Neoplasms From Ovarian Mucinous Neoplasms. Am J Clin Pathol. 2017 May 1;147(5):484-491.

6: Falconieri G, Cataldi P, Kavalari R, Štitič V, Luzar B. Cutaneous Osteoblastic Osteosarcoma: Report of 2

New Cases Integrated With SATB2 Immunohistochemistry and Review of the Literature. Am J Dermatopathol. 2016 Nov;38(11):824-831.

7: Li M, Cai YP, Lu KY, Chen Y, Zhu X, Yin Y, Tang J. [Significance of SATB2 in the pathologic diagnosis of osteosarcoma]. Zhonghua Bing Li Xue Za Zhi. 2016 Sep 8;45(9):631-5.

8: Machado I, Navarro S, Picci P, Llombart-Bosch A. The utility of SATB2 immunohistochemical expression in distinguishing between osteosarcomas and their malignant bone tumor mimickers, such as Ewing sarcomas and chondrosarcomas. Pathol Res Pract. 2016 Sep;212(9):811-6.





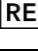
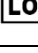
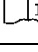
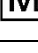

11: Sangoi AR, Kshirsagar M, Horvai AE, Roma AA. SATB2 Expression Is Sensitive but Not Specific for Osteosarcomatous Components of Gynecologic Tract Carcinosarcomas: A Clinicopathologic Study of 60 Cases. Int J Gynecol Pathol. 2017 Mar;36(2):140-145.

12: Fukuhara M, Agnarsdóttir M, Edqvist PH, Coter A, Ponten F. SATB2 is expressed in Merkel cell carcinoma. Arch Dermatol Res. 2016 Aug;308(6):449-54.

13: Kim CJ, Baruch-Oren T, Lin F, Fan XS, Yang XJ, Wang HL. Value of SATB2 immunostaining in the distinction between small intestinal and colorectal adenocarcinomas. J Clin Pathol. 2016 Dec;69(12):1046-1050.

LABEL AND BOX SYMBOLS

Explanation of the symbols of the product label and box:

| | |
|---|--|
|  | Expiration date |
|  | Temperature limit |
|  | Manufacturer |
|  | Sufficient content for <n> assays |
|  | Catalog number |
|  | Lot code |
|  | Refer to the instructions of use |
|  | Medical product for <i>in vitro</i> diagnosis. |
|  | Material safety data sheet |